



# Better Counts for Better Science

## LUNA FX7™ Automated Cell Counter







# Beyond Cell Counts

The LUNA-FX7 provides the flexibility to meet a multitude of cell counting needs. In addition to the need for accurate and reliable counts, bio-pharmaceutical and manufacturing facilities face unique challenges with regards to meeting QA/QC criteria, user management, and data security compliance. The LUNA-FX7's new features such as the Bioprocess Mode, Quality Control Mode, and the CountWire Software will help companies and facilities meet and surpass quality and regulatory requirements.


## Multiple Slide Options, Simple Handling

**Versatile Slide Options.** The LUNA-FX7™ can accommodate multiple slide options. The multi-channel pipette compatible 8-channel slides maximize throughput; the 3-channel slides facilitate easier data handling and organization; and the single channel slides, with a nearly 5 µL counting volume, promote increased accuracy and precision. With 'One and Done' slide insertion, a slide is inserted once and the LUNA-FX7™ will count all the chambers without the need to remove and reinsert the slide.

	NEW LUNA™ 1-Channel Slides	LUNA™ Cell Counting Slides / PhotonSlides™	NEW LUNA™ 3-Channel Slides	NEW LUNA™ 8-Channel Slides
Compatible slides				
Sample throughput	1 sample	Up to 2 samples	Up to 3 samples	Up to 8 samples
Sample loading volume	50 µL	10 µL / chamber	10 µL / chamber	10 µL / chamber
Maximum analysis volume	5.1 µL	1.3 µL / chamber	1.3 µL / chamber	0.5 µL / chamber

## Unique Bioprocess Mode

**Cell Monitoring.** The LUNA-FX7™ offers a Bioprocess mode that monitors cell health and growth and tracks bioprocessing activities. By easily tagging a batch to a specific protocol, the progression of multiple batches may be monitored. The Bioprocess application calculates doubling times and charts proliferation curves, allowing you to trace the production and health of your cells. With up to a terabyte of internal storage, the LUNA-FX7™ can store days, weeks, months, and years of counting results. Cell number concentrations and viabilities are provided with a single tap on the screen.

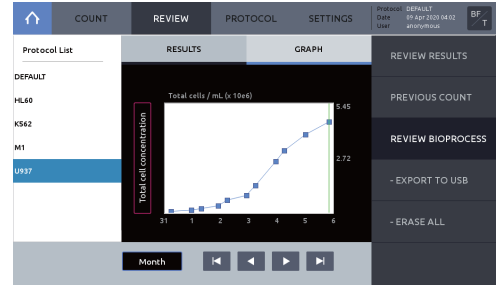

The screenshot shows the 'RESULTS' tab of the software interface. It displays detailed data for a specific protocol (L1937):

Parameter	Value
Start	31 Mar 2020 05:21:59
Runs	11
Last result	
Date	05 Apr 2020 21:11:13
Elapsed	135.821 hours
Growth rate	0.02189 / hour
Doubling time	22.7126 hours
Total cells	4.54 x 10 <sup>6</sup> cells / mL
Viability	24.1 %

## Quality Control - Trust in Your Data

**Innovation at Your Fingertips.** Our innovative quality control feature ensures that you can be confident in your data. In conjunction with our pre-fixed fluorescent or brightfield validation slides, the QC functions allow you to store, monitor, and analyze instrument performance over time. By establishing user-defined acceptance limits, the automated tracking and charting functions will help ensure compliance with the most rigorous QA/QC criteria.



# CountWire™-FDA 21 CFR Part 11

**FDA Regulation Ready.** Established by the U.S. Food & Drug Administration, 21 CFR Part 11, defines the criteria under which electronic records and electronic signatures are considered trustworthy, reliable, and equivalent to paper records. The FDA 21 CFR Part 11 regulations focus on the security and confidentiality of information, administration of e-documents, and the approval requirements for e-documents and signatures.


The LUNA-FX7™ is equipped with the innovative **CountWire™ software**. With user-accessibility management including encrypted activity logs, it delivers data securely and remains compliant with 21 CFR Part 11 regulations and GMP requirements.



Date	Device	User	File	User	Supervisor	Admin
20110024_193310	external	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	192.168.0.210@Bio.com	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	192.168.0.210@Bio.com	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	192.168.0.210@Partners 1	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	192.168.0.210@Partners 2	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	external	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	192.168.0.210@Bio.com	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	external	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	192.168.0.210@Partners 2	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	external	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	external	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	external	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	external	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	external	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	external	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	external	WJensen	210001-Std-001	WJensen	WJensen	WJensen
20110024_193310	external	WJensen	210001-Std-001	WJensen	WJensen	WJensen

Date	Device	User	File
20110024_193310	external	WJensen	210001-Std-001
20110024_193310	external	WJensen	210001-Std-001

**The CountWire™: Remote Access.** The CountWire™ software enables remote access and data management for multiple LUNA-FX7™s connected to your network. With the right level of user access, you can retrieve data from the LUNA-FX7™ in the clean room next door without leaving your desk.

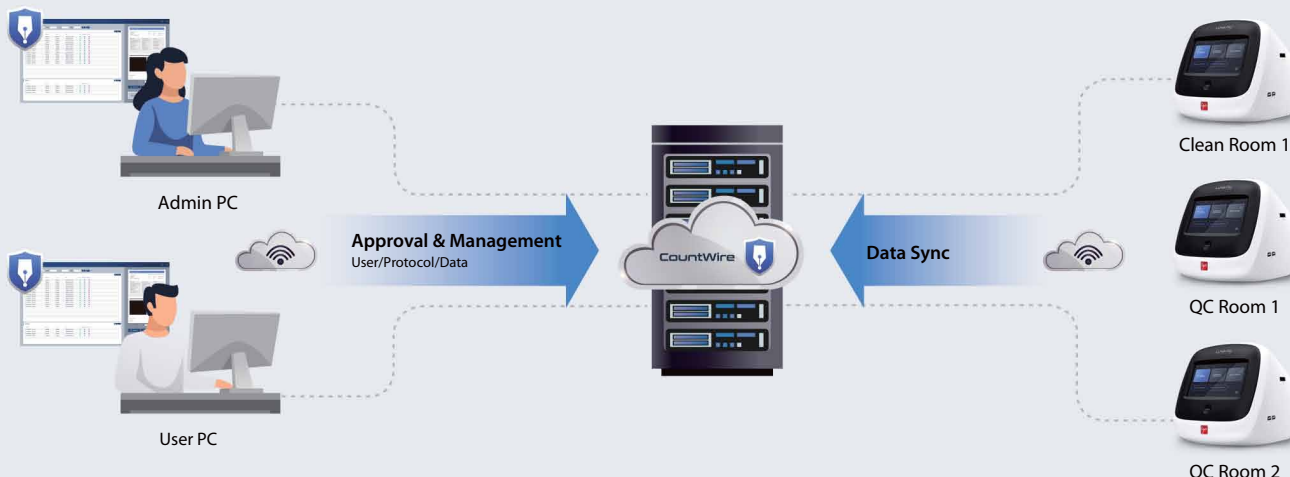

# Everything you love about our LUNAs and more

Unmatched cell counting accuracy

Optimized for bioprocess production applications

21 CFR Part 11 ready

- ✓ All-new optics
- ✓ Increased counting volume for the lowest CV per count
- ✓ Multichannel pipette-compatible 8-channel slides
- ✓ Fast and precise autofocus
- ✓ More robust and sophisticated counting algorithms
- ✓ Customizable cell-detection protocols
- ✓ Quality control and validation software
- ✓ Range of standard validation slides
- ✓ User access and rights management
- ✓ Online data storage and control
- ✓ Encrypted electronic records

## SPECIFICATIONS

LUNA-FX7™ Automated Cell Counter L70001



LUNA-FX7™ Automated Cell Counter Bioprocess L70002



Onboard storage	250 GB	1 TB
Software	Standard software	Bioprocess software
Cell diameter range	Detectable: 1-90 µm, optimal: 3-60 µm	
Optics	Brightfield, dual fluorescence	
Green fluorescence	Ex : 470±20 nm, Em : 530±25 nm	
Red fluorescence	Ex : 530±20 nm, Em : 620±30 nm	
Focusing	Auto focus / Manual focus	
Display	7" Capacitive touch screen LCD	
Data format	PDF, CSV, TIFF	
Data export	USB, WiFi, Ethernet	
Printer	External printer (optional)	
21 CFR Part 11	CountWire™ software (optional)	
User management	CountWire™ software (optional)	
IQ/OQ	Yes (optional)	Yes
Dimensions	245 x 280 x 240 mm	
Weight	4.5 kg	

## CELL LINES VALIDATED On The LUNA-FX7™ Automated Cell Counter

Cell Type	Vendor	Animal	Organ
3T3-L1	ATCC	Mouse	Embryo
A549	ATCC	Human	Lung
AsPC-1	ATCC	Human	Pancreas
BAEC	ATCC	Cow	Aorta
Colon26L5	ATCC	Mouse	Colon carcinoma
DLD-1	ATCC	Human	Colon
H1975	ATCC	Human	Lung
HCT116	ATCC	Human	Colon
Hek-293	ATCC	Human	Embryonic kidney
HL60	ATCC	Human	Peripheral blood
J774A.1	ATCC	Mouse	Ascites
RAW264.7	ATCC	Mouse	Abelson murine leukemia virus-induced tumor; ascites
SY5Y	ATCC	Human	Bone marrow
U2OS	ATCC	Human	Bone
U937	ATCC	Human	Pleura/pleural effusion, lymphocyte, Myeloid
K562	ATCC	Human	Bone marrow
PBMC	ATCC	Human	Blood
M1	ATCC	Mouse	Myeloblast

This is partial list of cell line tested on the LUNA-FX7™.

## ORDERING INFORMATION

LUNA-FX7™		L70001	LUNA-FX7™ Automated Cell Counter	1 unit	
		L70002	LUNA-FX7™ Automated Cell Counter, Bioprocess Package	1 unit	
Counting Slides	1-channel slides 50 µL <i>FX7-only</i>	L72011	LUNA™ 1-Channel Slides, 50 Slides	1 box	
		L72012	LUNA™ 1-Channel Slides, 500 Slides	10 boxes	
		L72013	LUNA™ 1-Channel Slides, Sterile - gamma-irradiated, 500 Slides	10 boxes	
	8-channel slides 10 µL <i>FX7-only</i>	L72001	LUNA™ 8-Channel Slides, 50 Slides	1 box	
		L72002	LUNA™ 8-Channel Slides, 500 Slides	10 boxes	
		L72003	LUNA™ 8-Channel Slides, Sterile - gamma-irradiated, 500 Slides	10 boxes	
	3-channel slides 10 µL <i>FX7-only</i>	L72021	LUNA™ 3-Channel Slides, 50 Slides	1 box	
		L72022	LUNA™ 3-Channel Slides, 500 Slides	10 boxes	
		L72023	LUNA™ 3-Channel Slides, Sterile - gamma-irradiated, 500 Slides	10 boxes	
	2-channel slides 10 µL <i>Brightfield</i>	L12001	LUNA™ Cell Counting Slides, 50 Slides	1 box	
		L12002	LUNA™ Cell Counting Slides, 500 Slides	10 boxes	
		L12003	LUNA™ Cell Counting Slides, 1000 Slides	20 boxes	
2-channel slides 10 µL <i>Fluorescence</i>	L12005	PhotonSlide™, 50 Slides	1 box		
	L12006	PhotonSlide™, 500 Slides	10 boxes		
	L12007	PhotonSlide™, 1000 Slides	20 boxes		
	Validation Slides		L72030	Cell Counter Validation Slide (Fluorescence)	1 unit
		L72040	Cell Counter Validation Slide (Brightfield)	1 unit	
Reagents			F73101	LUNA-FX7™ Calibration Beads Kit	1 kit
			T13001	Trypan Blue Stain, 0.4%	2 x 1 mL
			T13011	Trypan Blue Stain, 0.4%, Sterile-filtered	2 x 1 mL
			L13002	Erythrosin B Stain	2 x 1 mL
			F23001	Acridine Orange/Propidium Iodide Stain	2 x 0.5 mL
			F23011	Acridine Orange/Propidium Iodide Stain, Sterile-filtered	2 x 0.5 mL
			L71001	CountWire™ Basic	1 set
Software			L71002	CountWire™, Additional PC License	1 unit
	IQ/OQ		L74003	LUNA-FX7™ IQ/OQ Protocol	1 copy
Printer-related Items			P17001	Thermal Printer	1 unit
			P17002	Thermal Printer Refills, Paper and Ribbon	2 x 5 rolls



### HEADQUARTERS

FL 2 & 3  
28 Simindaero 327beon-gil, Dongan-gu  
Anyang-si, Gyeonggi-do 14055  
South Korea  
Email: [info@logosbio.com](mailto:info@logosbio.com)  
Tel: +82 (31) 478-4185  
Fax: +82 (31) 360-4277

### USA

7700 Little River Turnpike STE 207  
Annandale, VA 22003  
USA  
Email: [info-usa@logosbio.com](mailto:info-usa@logosbio.com)  
Tel: +1 (703) 622-4660, +1 (703) 942-8867  
Fax: +1 (571) 266-3925

### EUROPE

11B avenue de l'Harmonie  
59650 Villeneuve d'Ascq  
France  
Email: [info-france@logosbio.com](mailto:info-france@logosbio.com)  
Tel: +33 (0)3 74 09 44 35  
Fax: +33 (0)3 59 35 01 98